



Comprehensive Traceability for the Life of the Sample

**BlockTracker and SlideTracker:
Automated Archiving Solutions**

Protect and Organize Pathology Samples for the Future

The BlockTracker and SlideTracker Systems are designed to automate the filing and tracking of your physical slides and blocks for the life of the sample. The flexible design allows these units to fit into your current workflow, while automating the manual task of block and slide sorting and filing.

- Full traceability for the life of the sample
- Scalable solution for high throughput workflows
- Reduce archiving space requirements with high density storage
- Automate organization of slides and blocks into final storage configuration
- Improve resource utilization by automating manual archiving steps.

Location Based Archiving


The BlockTracker and SlideTracker, photographs and scans each sample and assigns it to a specific tray position. The magazine is then assigned to a storage location for traceability and easy retrieval.

Chain of Custody for Sample Traceability

Check in/out controls provides visibility to the user that retrieved the sample and why. Samples that are not returned in a defined timeframe will be flagged.

Automated Re-filing

Returned samples can be placed in any position and scanned into the system. The database will update the sample location automatically.

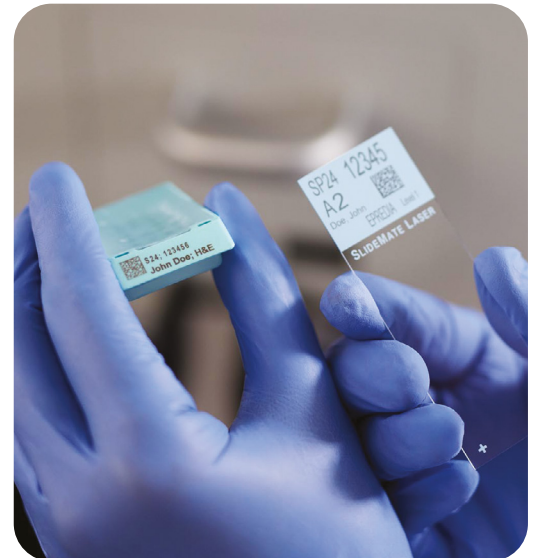

SlideTracker – Minimize hands on time for slide archiving

Scalable and automated workflow

- Archive up to 400 slides per hour automatically
- Eliminate manual filing of slides individually and the need to put in sequential order
- Automatically segments slides into different magazines based on the barcode
- Gain full traceability to every slides location and who entered the slide into the database
- Eliminate the need to refile slides into the same location, through automated updating of the database when a slide is returned.

Sample barcode allows automated identification of the location of each sample



SlideScan: On demand scanning for damaged or broken slides.



Slide magazine designed to maximize storage density and save up to 30% space.¹

¹ Data on file

Sample Management Software



Take control of your pathology assets

Consolidate block and slide management in a single automated software

Database software, allows users to manage all samples and gain full traceability to the movement and location of all blocks and slides that have been archived using Blocktracker and Slidetracker. Users can manage the location of slides and blocks from any PC connected to the network.

Samples can be automatically requested from offsite locations or from long-term storage. Consumable usage and future requirements are tracked and ordering patterns are suggested.

- Manage your sample archive in one user friendly software
- Quickly locate any sample in the database for retrieval or transfer of location
- Automate the picklist process through integration with your LIS or HIS system
- Access and manage the database from any PC on your network
- Run Reports on usage, storage location, consumable usage and more
- Backup your data locally or on the cloud to ensure your database is secure.

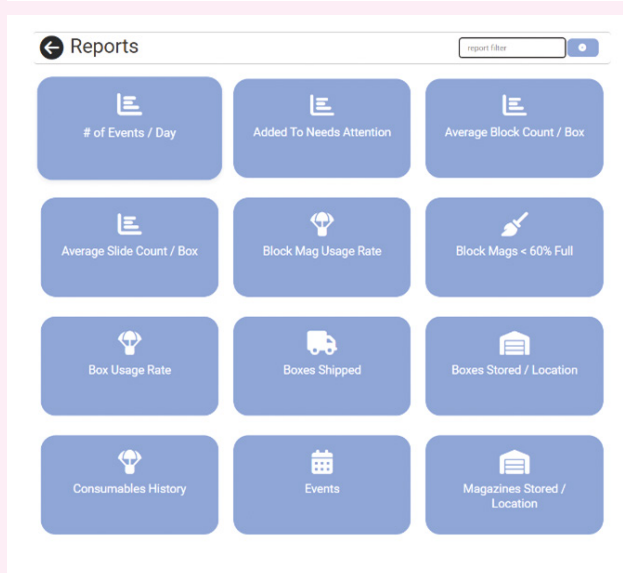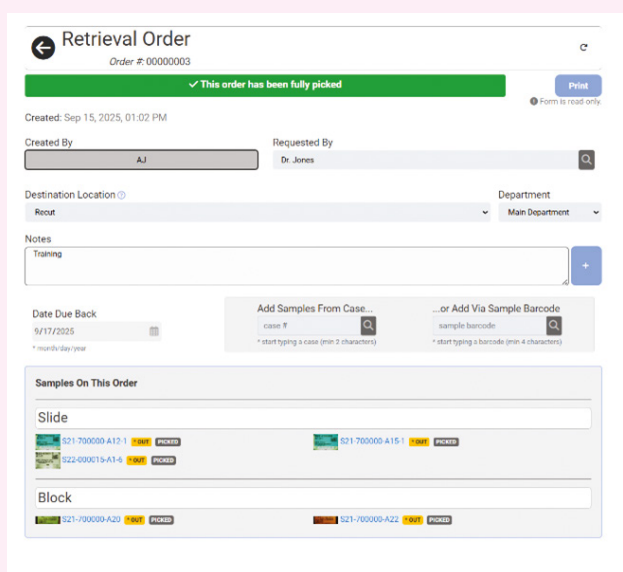



BlockTracker– Scalable solution for automating block archiving

- Scan and Archive up to 1500 blocks per hour
- Place blocks directly into storage magazines at the microtome and eliminate multiple manual filing steps
- Eliminate the need to place blocks in sequential and case order
- Automatically scan and photograph each block, to have full traceability of location and entry into the database
- Unreadable barcodes are segmented to alert the user with no disruption to the workflow
- Dense storage design allows you to maximize storage space and fits into current archiving configurations.





Specifications and Ordering Information

Cat. Number	Description	Capacity
STSYSTEM	SlideTracker System	400 slides/hour
BTSYSTEM	BlockTracker System	1500 blocks/hour
STBTSYSTEM	Combination BlockTracker and SlideTracker System	
SLIDESCAN	SlideScan Unit	
STMAG	SlideTracker Magazine Case QTY 96	200 slides per magazine
BTMAG	BlockTracker Magazine Case QTY 96	50 blocks per magazine

To find out about all laboratory solutions from Epredia, contact your Epredia Account Executive.
To place an order, contact Epredia Customer Care at ep-cs-kalamazoo@epredia.com; (800) 522-7270.

Find out more at [epredia.com](https://www.epredia.com)

4481 Campus Drive | Kalamazoo, MI 49008 | United States | +1 (800) 522-7270

© 2025 Epredia Group. All rights reserved. All trademarks are the property of Epredia Group and its subsidiaries unless otherwise specified. **M2500085 R072026**

