



Accuracy made easy

SlideMate AS Slide Printer




## How Does It Work?

- Data is sent from the LIS or scanned from the histology cassette to generate a barcode and slide identifier
- Templates can be customised
- Slides are printed using the intuitive software

## Slide Delivery System

- The optional Slide Delivery System is designed to hold 72 slides to reduce handling
- EpreDia SlideMate slides can be loaded using optimised packaging

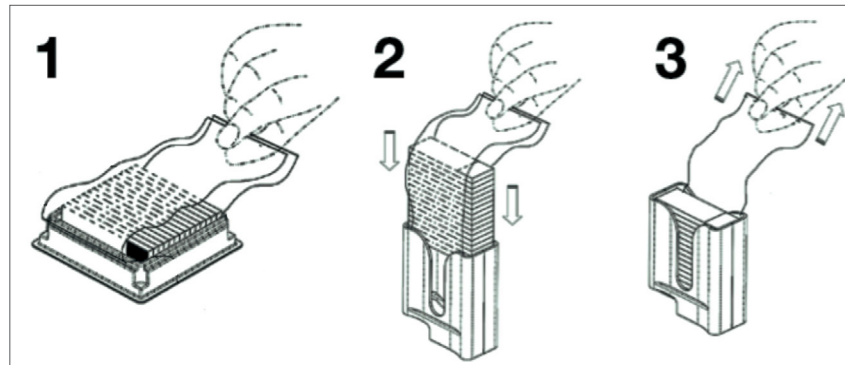

## Printing with SlideMate AS is easy as...



# SlideMate AS On-Demand Slide Printer

## Accurate Identification

- On-demand slide printing at the microtome can increase accuracy as patient identification is maintained

## Efficiency

- The SlideMate AS is designed to print slides within 5 seconds, enhancing workflow at the microtome

## Connectivity

- SlideMate AS has the option of being connected to other computers as well as LIS or LIMS using a Windows® driver

## Workflow Enhancement

- SlideMate AS may be used for barcoding slides within the workflow which allows for additional applications of automated archiving and digital pathology
- With its small footprint, the SlideMate AS slide printer fits directly next to your microtome

## The essential partner to your microtome



Fully automated HM355S rotary microtome

Optional slide hopper allows for quicker and easier on-demand printing of slides. Slides can automatically feed from the hopper.

Visit [www.epredia.com](http://www.epredia.com) to learn about SlideMate Pro dual hopper on-demand slide printer.



## SlideMate AS Printer Ordering Information

Systems	Order Number
<b>SlideMate AS Basic Printer</b> SlideMate AS basic printer with color touchscreen, power supply, power cord and Windows® Printer Driver software to print externally generated slide label from an image file	<b>A83900001</b>
<b>SlideMate AS Basic Printing System</b> SlideMate AS basic printer with slide delivery system, color touchscreen, power supply, power cord and Windows® Printer Driver software to print externally generated slide label from an image file	<b>A83900002</b>
<b>SlideMate AS On-Demand Printer</b> SlideMate AS on-demand printer with color touch-screen, built-in scanner, power supply, power cord and advanced software enabling label creation, special panel printing and import from LIS	<b>A83900004</b>
<b>SlideMate AS On-Demand Printing System</b> SlideMate AS on-demand printer with slide delivery system, color touch-screen, built-in scanner, power supply, power cord and advanced software enabling label creation, special panel printing and import from LIS	<b>A83900005</b>
Accessories	Cat. Number
Thermal transfer ribbon	A83910001
Slide delivery system	A83910002
Spare slide hopper for SlideMate AS slide delivery system	A83920016
Gateway software (1 unit)	A83910005
Multi-Gateway software (up to 10 units)	A83910009
Head cleaning kit	A83910008

## SlideMate AS Printer Specifications


Dimensions: SlideMate AS printer	
D x W X H	230 x 122 x 290 mm 9.1 x 4.8 x 11.5 in
Dimensions: SlideMate AS printer with slide delivery system and hopper	
D x W X H	230 x 122 x 470 mm 9.1 x 4.8 x 18.5 in
Printer Weight	4 kg, 9.5 lbs
Power requirements:	Universal 100 – 240 VAC, 50 – 60 Hz, 1.2A Max.
Slide type	25 mm x 75 mm x 1mm or 1 in x 3 in
Print speed	Up to 9 slides per minute with slide delivery system Up to 12 slides per minute with manual feed
Media type	Thermal transfer ribbon
Media capacity	1 roll, approximately 12,000 slides per roll
Media color	Black
Connectivity	Ethernet, USB-A (2), USB-B (1)

## SlideMate slides and SlideMate PLUS adhesion slides from EpreDia

Experience crisp, clean printing with SlideMate microscope slides from EpreDia, designed for the SlideMate AS slide printer.

SlideMate	Cat. Number
White	TT-40418217-W
Blue	TT-50418217-B
Green	TT-60418217-G
Pink	TT-70418217-P
Yellow	TT-80418217-Y

SlideMate PLUS	Cat. Number
White	TT-40418218-PS-W
Blue	TT-50418218-PS-B
Green	TT-60418218-PS-G
Pink	TT-70418218-PS-P
Yellow	TT-80418218-PS-Y

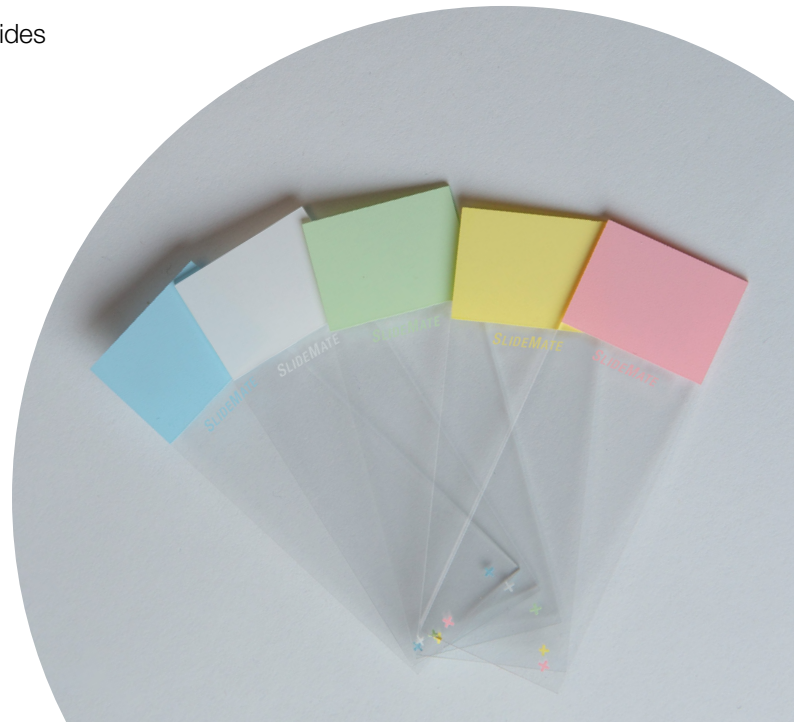

Find out more at [epredia.com](https://epredia.com)

4481 Campus Drive | Kalamazoo, MI 49008 | United States | +1 (800) 522-7270

© 2023 EpreDia Group. All rights reserved. All trademarks are the property of EpreDia and its subsidiaries unless otherwise specified. Windows® is a registered trademark of Microsoft®. These products are not available in all regions. Please reach out to your local sales representative for more information. **M41111\_EMEA\_A4\_R0623**

**epredia** Enhancing precision  
cancer diagnostics

POWERED BY SHANDON MICRON MENZEL-GLASER